

Data Sheet

LX-2 Human Hepatic Stellate Cell Line

Immortalized Cell Line

SCC064

Pack Size: 1x10⁶ viable cells/vial

Store in liquid nitrogen.

FOR RESEARCH USE ONLY

Not for use in diagnostic procedures. Not for human or animal consumption.

Background

Hepatic stellate cells are a major cell type responsible for liver fibrosis following their activation into fibrogenic myofibroblast-like cells in diseases such as chronic alcoholism, hepatitis B and C, fatty liver disease, obesity and diabetes. There is an increasing need for renewable cell culture models that faithfully recapitulate there in vivo phenotype, particularly for human studies.

LX-2 was generated by immortalization of primary human hepatic stellate cells with the SV40 large T antigen followed by selective culture of early passaged cells in low serum media conditions.

Immortalized LX-2 was established by Xu et al to overcome issues of culture variability and to provide a stable and unlimited source of human hepatic stellate cells that are homogeneous. These cell lines have been extensively characterized and retain key features of cytokine signaling, neuronal gene expression, retinoid metabolism, and fibrogenesis, making them highly suitable for culture-based studies of human hepatic fibrosis.

Source

Human hepatic stellate primary cells were isolated from a consenting normal human donor following established protocols outlined by Friedman et al.

Quality Control Testing


- Each vial contains $\geq 1 \times 10^6$ viable cells
- Cells are tested by PCR and are negative for Hepatitis A, B, C and HIV-1 & 2 viruses
- Cells are negative for mycoplasma contamination

Storage and Handling


LX-2 cells should be stored in liquid nitrogen. The cells can be passage for at least 10 passages without significantly affecting the cell marker expression and functionality. LX-2 cells have been successfully expanded past passage 50 in culture.

Representative Data


LX-2 Cells 24 hrs after thaw

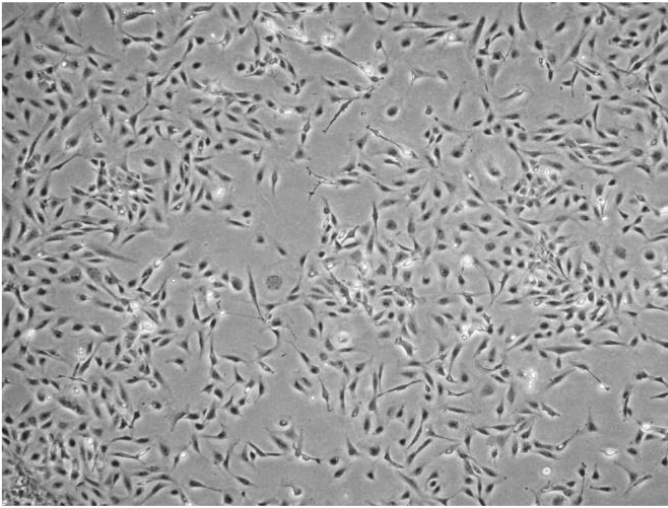


LX-2 Cells 72 hrs after thaw

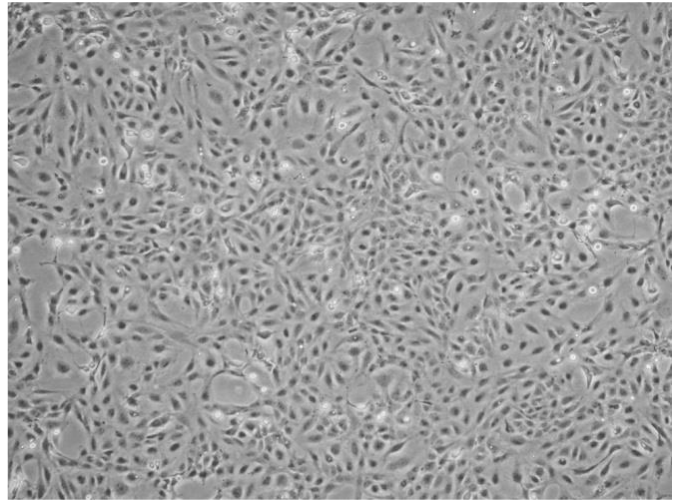


LX-2 Cells staining with α -SMA (1:50, MAB1501X)

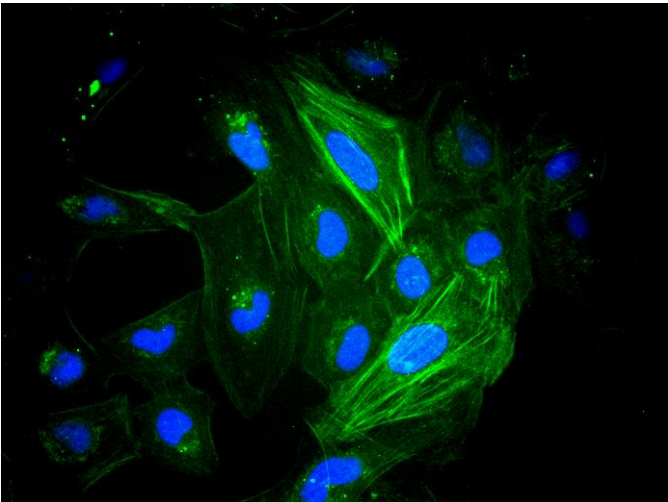

LX-2 Cells staining with Vimentin (10 μ g/mL, MAB3400)

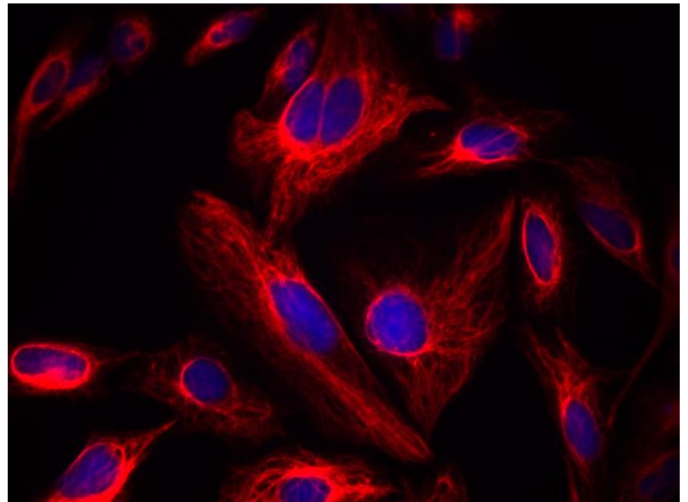

LX-2 Cells staining with GFAP (1:1000, AB5804)

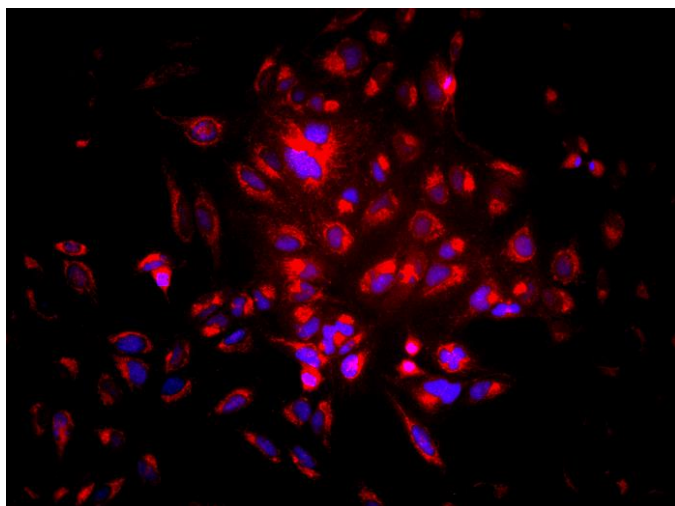

Protocols

Thawing of Cells

1. Do not thaw the cells until the recommended medium is on hand. Cells can grow on normal tissue culture ware surfaces without any additional coating. Cells are thawed in DMEM High Glucose (D5796), 10% FBS (ES009-B), 1X Pen/Strep (TMS-AB2-C) and 1X Glutamine (TMS-002-C) media. Once thawed, cells are expanded in 2% FBS media using the same components listed above.
2. Remove the vial of LX-2 cells from liquid nitrogen and incubate in a 37 °C water bath. Closely monitor until the cells are completely thawed. Maximum cell viability is dependent on the rapid and complete thawing of frozen cells.
IMPORTANT: Do not vortex the cells.
3. As soon as the cells are completely thawed, disinfect the outside of the vial with 70% ethanol. Proceed immediately to the next step.
4. In a laminar flow hood, use a 1- or 2- mL pipette to transfer the cells to a sterile 15 mL conical tube. Be careful not to introduce any bubbles during the transfer process.
5. Using a 10 mL pipette, slowly add dropwise 9 mL of 10% FBS media (Step 1 above) (pre-warmed to 37 °C) to the 15 mL conical tube.
IMPORTANT: Do not add the whole volume of media at once to the cells. This may result in decreased cell viability due to osmotic shock.
6. Gently mix the cell suspension by slow pipeting up and down twice. Be careful to not introduce any bubbles.
IMPORTANT: Do not vortex the cells.
7. Centrifuge the tube at 300 x *g* for 2-3 minutes to pellet the cells.
8. Decant as much of the supernatant as possible. Steps 5-8 are necessary to remove residual cryopreservative (DMSO).
9. Resuspend the cells in a total volume of 10 % FBS medium (prewarmed to 37 °C).
10. Plate the cell mixture onto a T75 tissue culture flask.
11. Incubate the cells at 37 °C in a 5% CO₂ humidified incubator.
12. The next day, exchange the medium with fresh 2% FBS media (Step 1 above) pre-warmed to 37 °C. Exchange with fresh medium every 2-3 days thereafter.
13. When the cells are approximately 80% confluent (3-4 days after plating cells at the density they can be dissociated with Accutase® (SCR005) or trypsin (SM-2003-C) and passaged or alternatively frozen for later use.

Subculturing of Cells

1. Carefully remove the medium from the T75 tissue culture flask containing the confluent layer of LX-2 cells.
2. Apply 3-5 mL of Accutase® or trypsin solution and incubate in a 37 °C incubator for 3-5 minutes.
3. Inspect the plate and ensure the complete detachment of cells by gently tapping the side of the plate with the palm of your hand.
4. Add 8 mL of 2% FBS medium (pre-warmed to 37 °C) to the plate.
5. Gently rotate the plate to mix the cell suspension. Transfer the dissociated cells to a 15 mL conical tube.
6. Centrifuge the tube at 300 x *g* for 3-5 minutes to pellet the cells.
7. Discard the supernatant.
8. Apply 2 mL of 2% FBS media (pre-warmed to 37 °C) to the conical tube and resuspend the cells thoroughly.
IMPORTANT: Do not vortex the cells.
9. Count the number of cells using a hemocytometer.
10. Plate the cells to the desired density (typical split ratio is 1:3 to 1:6).

Cryopreservation of Cells

LX-2 cells grown in 2% FBS media can be frozen in 20% FBS growth media plus 10% DMSO using a Nalgene® slow freeze Mr. Frosty® container.

References

1. Xu, L., et al. (2005) Human hepatic stellate cell lines, LX-1 and LX-2: new tools for analysis of hepatic fibrosis. *Gut* 54 (1):142–51.
2. Friedman, S. L., et al. (1992) Isolated hepatic lipocytes and Kupffer cells from normal human liver: morphological and functional characteristics in primary culture. *Hepatology* 15(2): 234-43.
3. Weiskirchen, R., et al. (2013) Genetic Characteristics of the Human Hepatic Stellate Cell Line LX-2. *PLoS One* 8(10):e75692.
4. Chen, L., et al. (2013) Epigenetic regulation of connective tissue growth factor by microRNA-214 delivery in exosomes from mouse or human hepatic stellate cells. *Hepatology* [Epub ahead of print].
5. Arriola Benitez, P.C., et al. (2013) *Brucella abortus* induces collagen deposition and MMP-9 down-modulation in hepatic stellate cells via TGF-β1 production. *Am J Pathol* 183(6):1918-27.
6. Feng, G., et al. (2013) Interleukin-21 mediates hepatitis B virus-associated liver cirrhosis by activating hepatic stellate cells. *Hepatol Res.* [Epub ahead of print].
7. Liao, R., et al. (2013) Clinical significance and gene expression study of human hepatic stellate cells in HBV related hepatocellular carcinoma. *J Exp Clin Cancer Res.* 32:22.
8. Park, S.J., et al. (2013) TRAIL regulates collagen production through HSF1-dependent Hsp47 expression in activated hepatic stellate cells. *Cell Signal* 25(7):1635-43.
9. Sun, X., et al. (2013) Expression of Septin4 in human hepatic stellate cells LX-2 stimulated by LPS. *Inflammation* 36(3):539-48.
10. Wu, C. G., et al. (2013) Hepatitis C virus core protein stimulates fibrogenesis in hepatic stellate cells involving the obese receptor. *J Cell Biochem* 114(3):541-50.

Academic Use Agreement

Subject to local law

THIS PRODUCT MAY ONLY BE USED BY INDIVIDUALS EMPLOYED BY AN ACADEMIC INSTITUTION AND IS INTENDED SOLELY TO BE USED FOR ACADEMIC RESEARCH, WHICH IS FURTHER DEFINED BELOW. BY OPENING THIS PRODUCT, YOU ("PURCHASER") HEREBY REPRESENT THAT YOU HAVE THE RIGHT AND AUTHORITY TO LEGALLY BIND YOURSELF AND/OR YOUR EMPLOYER INSTITUTION, AS APPLICABLE, AND CONSENT TO BE LEGALLY BOUND BY THE TERMS OF THIS ACADEMIC USE AGREEMENT. IF YOU DO NOT AGREE TO COMPLY WITH THESE TERMS, YOU MAY NOT OPEN OR USE THE PRODUCT AND YOU MUST CALL CUSTOMER SERVICE (1-800-645-5476) TO ARRANGE TO RETURN THE PRODUCT FOR A REFUND.

"Product" means LX-2 Human Hepatic Stellate Cell Line (SCC064).

"Academic Research" means any internal in vitro research use by individuals employed by an academic institution. Academic Research specifically excludes the following uses of whatever kind or nature:

- Re-engineering or copying the Product
- Making derivatives, modifications, or functional equivalents of the Product
- Obtaining patents or other intellectual property rights claiming use of the Product
- Using the Product in the development, testing, or manufacture of a Commercial Product
- Using the Product as a component of a Commercial Product
- Reselling or licensing the Product
- Using the Product in clinical or therapeutic applications including producing materials for clinical trials
- Administering the Product to humans
- Using the Product in collaboration with a commercial or non-academic entity

"Commercial Product" means any product intended for: (i) current or future sale; (ii) use in a fee-for-service; or (iii) any diagnostic, clinical, or therapeutic use.

Access to the Product is limited solely to those officers, employees, and students of PURCHASER's academic institution who need access to the Product to perform Academic Research. PURCHASER shall comply with all applicable laws in its use and handling of the Product and shall keep it under reasonably safe and secure conditions to prevent unauthorized use or access.

These use restrictions will remain in effect for as long as PURCHASER possesses the Product.

COMMERCIAL OR NON-ACADEMIC ENTITIES INTERESTED IN PURCHASING OR USING THE PRODUCT MUST CONTACT licensing@emdmillipore.com AND AGREE TO SEPARATE TERMS OF USE PRIOR TO USE OR PURCHASE.

Genetically Modified Organisms (GMO)

This product contains genetically modified organisms.
Este producto contiene organismos genéticamente modificados.
Questo prodotto contiene degli organismi geneticamente modificati.
Dieses Produkt enthält genetisch modifizierte Organismen.
Ce produit contient des organismes génétiquement modifiés.
Dit product bevat genetisch gewijzigde organismen.
Tämä tuote sisältää geneettisesti muutettuja organismeja.
Denna produkt innehåller genetiskt ändrade organismer.

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck KGaA, Darmstadt, Germany
operates as MilliporeSigma in the U.S. and Canada.

MilliporeSigma, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2008-2025 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

Document Template 20306518 Ver 6.0

SCC064CA Ver 4.0, Rev 01JAN2025, AV

**Millipore
Sigma**