


Product Information

D-Amphetamine hemisulfate salt solution
1.0 mg/mL $\pm$ 5% in methanol
analytical standard, for drug analysis

Catalog Number **A3278**
Lot Number 041M8735
Storage Temperature 2–8 °C

CAS RN 51-63-8
Synonyms: Dextroamphetamine hemisulfate salt;
(+)- α -Methylphenethylamine hemisulfate salt;
(+)- α -Methylbenzeneethanamine hemisulfate;
Phenyl-2-aminopropane hemisulfate; Phenaminum

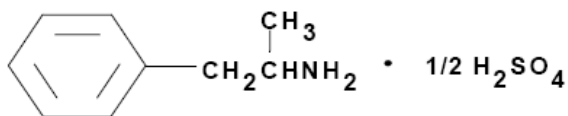

Product Description

Molecular Formula: $C_9H_{13}N \cdot 0.5 H_2SO_4$
Molecular Weight: 184.25

DEA Class: Exempt preparation of Class II

Concentration: 1.35 mg of D-Amphetamine hemisulfate (equivalent to 0.99 mg of free base)/mL of GC grade methanol, concentration verified by UV.

Expiration Date: July 2014

Precautions and Disclaimer


This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

Store the product at 2–8 °C. If stored at 2–8 °C in an airtight container, protected from light, decomposition of the product is less than 1% in 48 months. After opening, concentration may change due to loss of solvent.

CMW,MAM 07/11-2

70 eV Electron Ionization Mass Spectrum

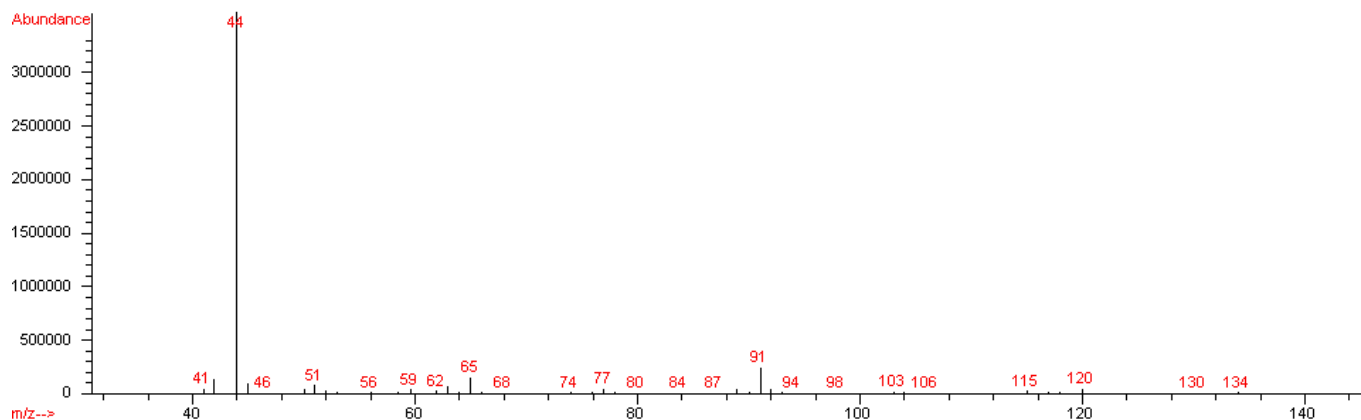

Hewlett-Packard 5890 Series II Gas Chromatograph
Hewlett-Packard 5971 Mass Selective Detector
Hewlett-Packard G1034C MS ChemStation (HP-DOS series)


FT-IR



Nicolet 380 FT-IR Thermo Smart Orbit

Ultraviolet Spectrum


Peak	258 nm
Absorbance	0.6604
Solvent	Methanol
Dilution	1:1 (v/v)

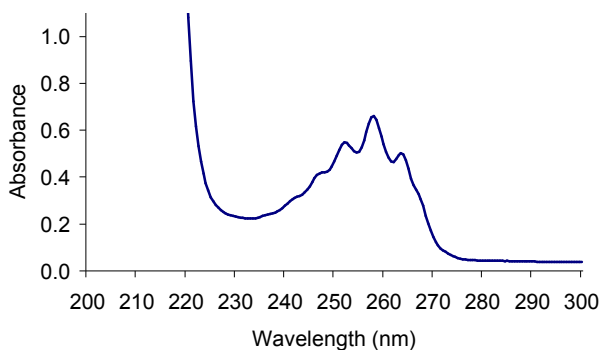

Beckman-Coulter DU 800

Capillary Gas Chromatography

Column Supelco SPB-1 (Catalog Number 24162)
I.D. 0.20 mm **Length** 15 m **Film Thickness** 0.20 μ m
Oven Temperature 50 °C (2 minutes), then
 15 °C/minute to 300 °C, hold 11 minutes
Injector Temperature 250 °C, Splitless
Detector Hewlett-Packard 5971
 Mass Selective Detector (70 eV, 40-500 m/z)

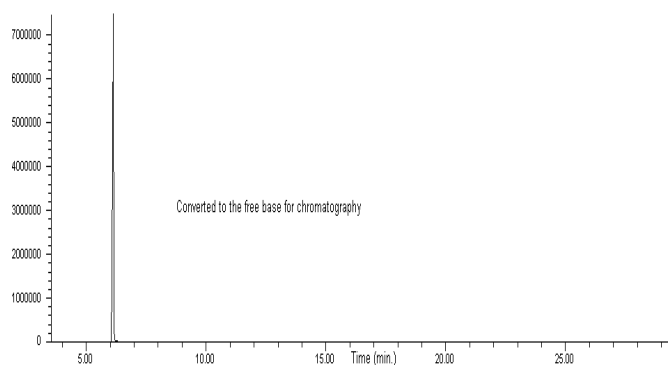

Hewlett-Packard 5890

Sigma brand products are sold through Sigma-Aldrich, Inc.

Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.