


## Product Information

### SIGMA QUALITY CONTROL TEST PROCEDURE

#### Enzymatic Assay of ENOLASE (EC 4.2.1.11)

##### PRINCIPLE:

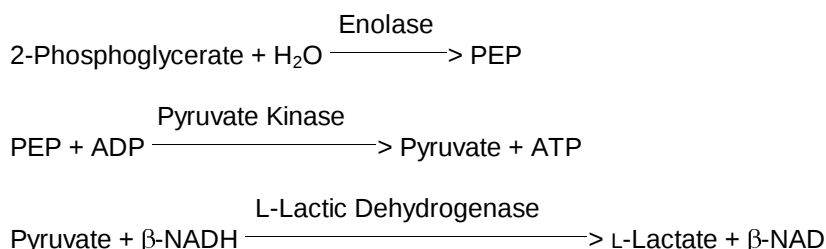

Abbreviations used:

PEP = Phospho(enol)Pyruvate

ADP = Adenosine 5'-Diphosphate

ATP = Adenosine 5'-Triphosphate

$\beta$ -NADH =  $\beta$ -Nicotinamide Adenine Dinucleotide, Reduced Form

$\beta$ -NAD =  $\beta$ -Nicotinamide Adenine Dinucleotide, Oxidized

**CONDITIONS:** T = 25°C, pH = 7.4,  $A_{340\text{nm}}$ , Light path = 1 cm

**METHOD:** Continuous Spectrophotometric Rate Determination

##### REAGENTS:

- A. 100 mM Triethanolamine Buffer, pH 7.4 at 25°C  
(Prepare 50 ml in deionized water using Triethanolamine Hydrochloride, Sigma Prod. No. T-1502. Adjust to pH 7.4 at 25°C with 1 M NaOH.)
- B. 56 mM 2-Phosphoglycerate Solution (DPG)  
(Prepare 2 ml in deionized water using D(+)-2-Phosphoglyceric Acid, Sodium Salt, Hydrate, Sigma Prod. No. P-0257.)
- C. 7 mM  $\beta$ -Nicotinamide Adenine Dinucleotide, Reduced Form Solution ( $\beta$ -NADH)  
(Dissolve the contents of one 5 mg vial of  $\beta$ -Nicotinamide Adenine Dinucleotide, Reduced Form, Disodium Salt, Sigma Stock No. 340-105, in 1 ml of Reagent A or prepare 1 ml in deionized water using  $\beta$ -Nicotinamide Adenine Dinucleotide, Reduced Form, Disodium salt, Sigma Prod. No. N-8129. **PREPARE FRESH.**)
- D. 500 mM Magnesium Sulfate with 2 M Potassium Chloride Solution ( $\text{MgSO}_4/\text{KCl}$ )  
(Prepare 5 ml in deionized water using Magnesium Sulfate, Anhydrous, Sigma Prod. No. M-7506, and Potassium Chloride, Sigma Prod. No. P-4504.)

## Enzymatic Assay of ENOLASE (EC 4.2.1.11)

### REAGENTS: (continued)

- E. 20 mM Adenosine 5'-Diphosphate Solution (ADP)  
(Prepare 1 ml in deionized water using Adenosine 5'-Diphosphate, Sodium Salt, Sigma Prod. No. A-2754. **PREPARE FRESH.**)
- F. PK/LDH Mixed Enzymes<sup>1</sup>  
(Use PK/LDH Enzyme Solution, Sigma Prod. No. P-0294)
- G. 15 mM Tris HCl with 0.02% (w/v) Bovine Serum Albumin (Enzyme Diluent)  
(Prepare 100 ml in deionized water using Trizma Base, Sigma Prod. No. T-1503, and Albumin, Bovine, A-4503, or equivalent. Adjust to pH 7.4 at 25°C with 1 M HCl.)
- H. Enolase Enzyme Solution  
(Immediately before use, prepare a solution containing 0.25 - 0.5 unit/ml of Enolase in cold Reagent G.)

### PROCEDURE:

Pipette (in milliliters) the following reagents into suitable cuvettes:

|                                          | <u>Test</u> | <u>Blank</u> |
|------------------------------------------|-------------|--------------|
| Reagent A (Buffer)                       | 2.39        | 2.39         |
| Reagent B (DPG)                          | 0.10        | 0.10         |
| Reagent C ( $\beta$ -NADH)               | 0.05        | 0.05         |
| Reagent D ( $\text{MgSO}_4/\text{KCl}$ ) | 0.15        | 0.15         |
| Reagent E (ADP)                          | 0.20        | 0.20         |
| Reagent F (PK/LDH)                       | 0.01        | 0.01         |

Mix by inversion and equilibrate to 25°C. Monitor the  $A_{340\text{nm}}$  until constant, using a suitably thermostatted spectrophotometer. Then add:

|                             |       |       |
|-----------------------------|-------|-------|
| Reagent G (Enzyme Diluent)  | ----- | 0.10  |
| Reagent H (Enzyme Solution) | 0.10  | ----- |

## Enzymatic Assay of ENOLASE (EC 4.2.1.11)

### PROCEDURE:(continued)

Immediately mix by inversion and record the decrease in  $A_{340\text{nm}}$  for approximately 5 minutes. Obtain the  $\Delta A_{340\text{nm}}$ /minute using the maximum linear rate for both the Test and Blank.

### CALCULATIONS:

$$\text{Units/ml enzyme} = \frac{(\Delta A_{340\text{nm}}/\text{min Test} - \Delta A_{340\text{nm}}/\text{min Blank})(3)(\text{df})}{(6.22)(0.1)}$$

3 = Volume (in milliliters) of assay

df = Dilution factor

6.22 = Millimolar extinction coefficient of  $\beta$ -NADH at 340 nm

0.1 = Volume (in milliliter) of enzyme used

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

### UNIT DEFINITION:

One unit will convert 1.0  $\mu\text{mole}$  of 2-phosphoglycerate to phospho(enol)pyruvate per minute at pH 7.4 at 25°C.

### FINAL ASSAY CONCENTRATION:

In a 3.00 ml reaction mix, the final concentrations are 81 mM triethanolamine, 1.9 mM 2-phosphoglycerate, 0.12 mM  $\beta$ -nicotinamide adenine dinucleotide, reduced form, 25 mM magnesium sulfate, 100 mM potassium chloride, 1.3 mM adenosine 5'-diphosphate, 7 units pyruvate kinase, 10 units L-lactic dehydrogenase and 0.025 - 0.05 unit Enolase.

### REFERENCES:

Bergmeyer, H.U. (1974) *Methods in Enzymatic Analysis*, 2nd ed., Vol. II, 449.

### NOTES:

1. Contains approximately 700 units/ml of Pyruvate Kinase and 1000 units/ml of Lactic Dehydrogenase.

**Enzymatic Assay of ENOLASE  
(EC 4.2.1.11)**

**NOTES:** (continued)

2. Pyruvate Kinase unit definition: One unit will convert 1.0  $\mu$ mole of phospho(enol)pyruvate to pyruvate per minute at pH 7.6 at 37°C.
3. Lactic Dehydrogenase unit definition: One unit will reduce 1.0  $\mu$ mole of pyruvate to L-lactate per minute at pH 7.5 at 37°C.
4. All product and stock numbers, unless otherwise indicated, are Sigma product and stock numbers.

**Sigma warrants that the above procedure information is currently utilized at Sigma and that Sigma products conform to the information in Sigma publications. Purchaser must determine the suitability of the information and products for its particular use. Upon purchase of Sigma products, see reverse side of invoice or packing slip for additional terms and conditions of sale.**