

Product Information

Tyrosinase

From mushroom

T3824

Product Description

CAS Registry Number: 9002-10-2

EC (Enzyme Commission) Number: 1.14.18.1

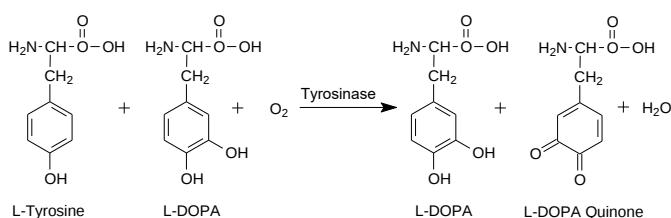
Synonyms: Monophenol monooxygenase; Polyphenol oxidase; Monophenol, dihydroxyphenylalanine:oxygen oxidoreductase;

Molecular mass:¹ 119.5 kDa (electrophoresis)Isoelectric point (pI):² 4.7–5

pH optimum: 6–7

Tyrosinase is a copper-containing oxidase which has activity on both catechols and cresol. Tyrosinase is responsible for browning reactions. Tyrosinase is reported to have two binding sites for aromatic substrates and a different binding site for oxygen, the copper-containing site.³

Overall Reaction

 K_M for L-Tyrosine:⁴ 0.5 mM K_M for L-DOPA:⁴ 1.5 mM

Catechol acts as a suicide substrate for tyrosinase, causing irreversible inactivation during catechol oxidation.⁵ Tyrosinase is also inhibited by compounds that complex with copper, such as benzoic acid and cyanide. Benzoic acid inhibition is competitive with catechol. Cyanide inhibition is competitive with oxygen and noncompetitive with catechol.³

This product is isolated from the mushroom species *Agaricus bisporus*.⁶ The product was dialyzed against deionized water prior to lyophilization. This enzyme is also assayed for polyphenol oxidase and catechol oxidase activities.

Activity: $\geq 1,000$ units/mg solid (tyrosinase activity)

Unit Definition: One unit will cause an increase in A_{280} of 0.001 per minute at pH 6.5 at 25 °C, in a 3 mL reaction mix containing L-tyrosine.

Several publications,⁷ theses⁸⁻¹² and dissertations¹³⁻¹⁷ have cited use of product T3824 in their research.

Precautions and Disclaimer

For R&D use only. Not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Preparation Instructions

Tyrosinase is soluble (2 mg/mL) in 50 mM potassium phosphate buffer, pH 6.5, yielding a clear brown solution.

Storage/Stability

When stored at -20 °C, the product loses <5% activity per year. If stored at 37 °C, 40-60% activity is lost in 7 days.¹⁸

Solutions retain activity for several days at 4-8 °C and for several weeks frozen at -20 °C.⁶ The enzyme is rapidly denatured in 1% SDS solution.¹ Tyrosinase dissociates into subunits in saturated urea and is completely inactivated. While the dissociation is reversible, the inactivation is not.⁶

References

1. Gillespie, J.P. *et al.*, *Comp. Biochem. Physiol. Pt. C*, **98(2-3)**, 351-358 (1991).
2. Robb, D.A., and Gutteridge, S., *Phytochemistry*, **20(7)**, 1481-1485 (1981).
3. Duckworth, H.W., and Colman, J.E., *J. Biol. Chem.*, **245(7)**, 1613-1625 (1970).
4. Espin, J.C. *et al.*, *J. Agric. Food Chem.*, **47(9)**, 3495-3502 (1999).
5. Cánovas, F.G. *et al.*, *Biochim. Biophys. Acta*, **912(3)**, 417-423 (1987).
6. Kertesz, D., and Zito, R., *Biochim. Biophys. Acta*, **96(3)**, 447-462 (1965).
7. Matsui, K.N. *et al.*, *J. Food Eng.*, **88(2)**, 169-176 (2008).
8. de la Torre, Jerry James M., "Effect of Frequency on Polyphenoloxidase Activity During Moderate Electric Field Treatment". The Ohio State University, M.S. thesis, pp. 22, 23 (2009).
9. Thallinger, Barbara, "Enzymatic crosslinking of human serum albumin with mouse IgG1 monoclonal antibody". Technischen Universität Graz, Dipl.-Ing. (M. Eng.) thesis, p. 15 (2011).
10. Thampi, Arya, "Unraveling a Few Mysteries of Melanin Polymerization: Structure and Kinetics of Formation". Indian Institute of Science Education and Research Pune, B.S./M.S. thesis, p. 6 (2015).
11. Adewinogo, Selena Eunice Orango, "Chemical Composition of the Essential Oils of Selected Aromatic South African Plants and Assessment of their Biological Activities for Cosmetic Uses". Cape Peninsula University of Technology, M.Sc. thesis, p. 53 (2020).
12. Prawitz, Lisa, "How Substances in Agricultural Waste Material Can Be Used as Active Ingredients in Skin Care Products". Lund University, M.Sc. thesis, p. 62 (2020).
13. Miyawaki, Miyuki, "Control of Polyphenol Oxidase and Pectin Methylesterase Activities by Ultra High Pressure". Washington State University, Ph.D. dissertation, p. 83 (2006).
14. Ismaya, Wangsa Tirta, "Crystal structure of *Agaricus bisporus* tyrosinase". University of Groningen, Ph.D. dissertation, p. 45 (2011).
15. Sarıçay, Yunus, "New Insight into Enzymatic Cross-linking of Globular Proteins: From Nanostructure to Functionality". Wageningen University, Ph.D. dissertation, p. 174 (2014).
16. Zakaria, Nik Nur Azwanida Binti, "Evaluation of anti-ageing properties of *Moringa oliefera* Lam., *Centella asiatica* (L.) Urban, *Clitoria ternatea* L. and *Cosmos caudatus* Kunth. for potential application as cosmeceuticals". Newcastle University, Ph.D. dissertation, p. 110 (2019).
17. Wallace, Andrea Kimi, "Engineering Diatom Peptides for the Synthesis of Silica Nanomaterials". Massachusetts Institute of Technology, Ph.D. dissertation, p. 105 (2020).
18. Kertesz, D. *et al.*, *Biochem. Biophys. Res. Commun.*, **49(5)**, 1208-1215 (1972).

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Standard Warranty

The applicable warranty for the products listed in this publication may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck operates as MilliporeSigma in the U.S. and Canada.

Merck and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

T3824dat Rev 03/22 RBG,GCY,MAM

